

## CLEANTYPE<sup>®</sup> PRIME SERIES

PANEL/PANEL+ ( WITH WALL MOUNT )



### SPECIAL FEATURES

- Foldable silicone control panel
- Adjustable brightness of backlight (PANEL+)
- 109 keys incl.num pad and track pad
- Lock after panel folded
- Excellent tactile feeling
- Integrated Clean Mode for cleaning
- EN60601-1-2 Ed 4 certified
- 2 USB extension port

### SPECIFICATION


Number of keys:	109(US-Layout)	Interfaces:	USB 2.0
Keyswitch technology:	carbon contact	Cable length:	150cm
Switching cycles:	approx.2 mill.(per key)	Operating temperature:	0°C~55°C
protection level:	IP65	Storage temperature:	-30°C~+70°C
Housing material:	silicone/stainless steel	Dimensions:	655(L)x215(W)x40(T)mm
Housing color:	metallic grey	Weight:	approx.3.9kg
Backlight:	blue color,5 level adjustable brightness	Certification:	EN60601-1-2 Ed 4,CE,FCC,RoHS
		Warranty:	2 years

### ORDER INFORMATION

Cat. No.	Product description	Color(keyboard)	Layout
PAL-U10220	CLEANTYPE <sup>®</sup> PRIME PANEL+	White	US/DE/JP/KR/THAI
PAL-U10221	CLEANTYPE <sup>®</sup> PRIME PANEL	White	US/DE/JP/KR/THAI
PAL-U10220	CLEANTYPE <sup>®</sup> PRIME PANEL+	Black	US/DE/JP/KR/THAI
PAL-U10221	CLEANTYPE <sup>®</sup> PRIME PANEL	Black	US/DE/JP/KR/THAI

Other layouts, configurations and interfaces on request.  
 Customized key function on request (Firmware update via USB)

## TECHNICAL DRAWING (IN MM)

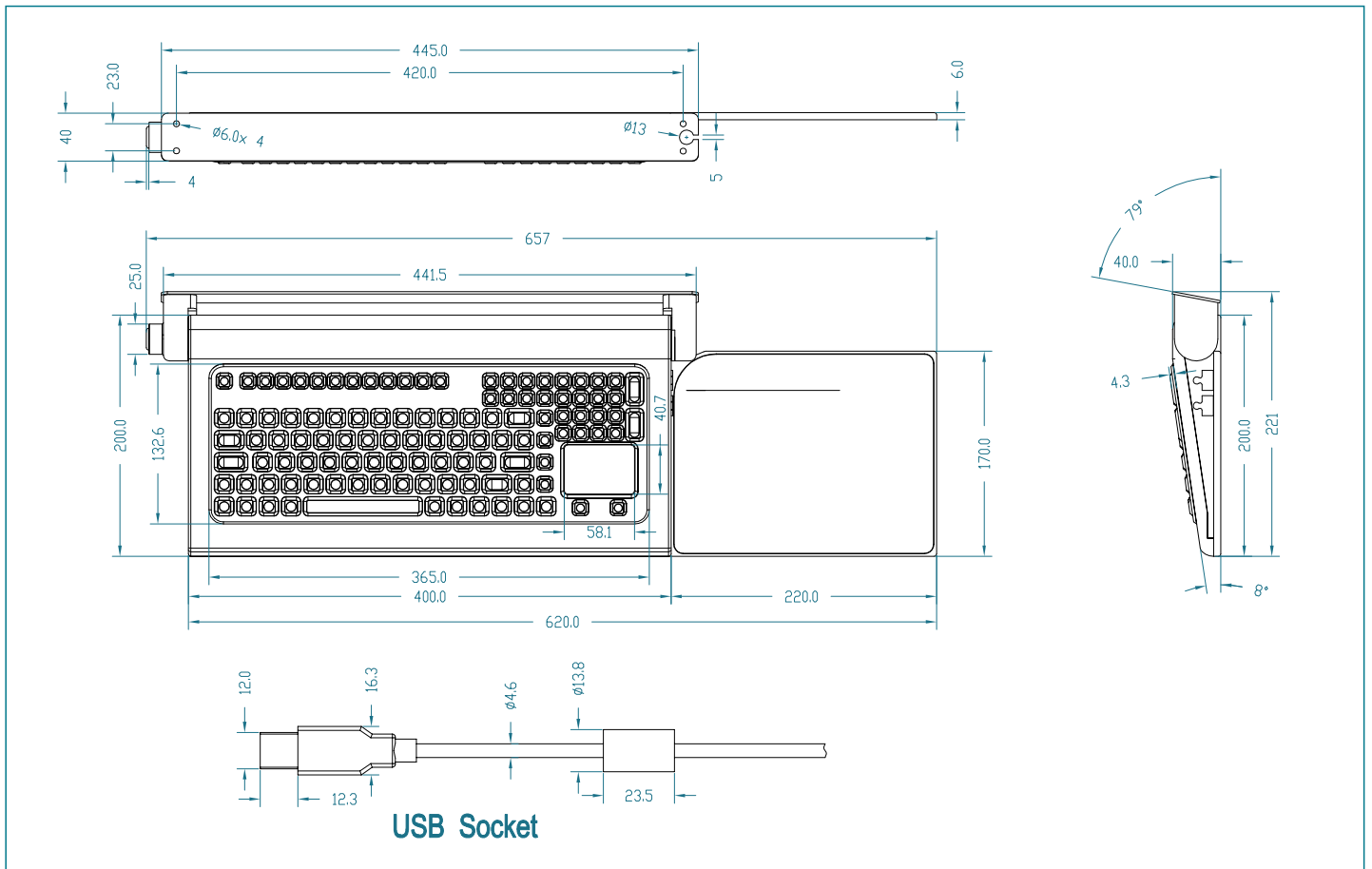

### EN60601-1-2 Ed 4 certified

EN 60601-1-2 Ed 4 applies to the basic safety and essential performance of Medical Equipment (ME) equipment and ME systems in the presence of electromagnetic disturbances and to electromagnetic disturbances emitted by ME equipment and ME systems.

### SHORT OPERATING INSTRUCTIONS



Press and hold lock key for 5s to enable Clean Mode. Clean Mode is indicated by breathing backlight.



Adjust backlight level with this button. Five level can be adjusted including off.

### CLEANING INSTRUCTION

For further details please refer to below link:  
<http://www.gett-group.com/content/service-cleaning-instructions>